

Extender™

The Extender™ is designed to resolve patient comfort issues that are typical with TLSOs trimmed above the breasts. Such TLSOs never seem to be set at the correct extension angle. By using HOPE Orthopedic's Extender to provide a thoracic extensive force, you achieve the necessary flexion control while enhancing patient compliance.

Features & Variations

- Sliding adjustable aluminum bars allow a custom fit
 - Easily attached to the pre-drilled front insert on the following HOPE spinal braces: ProLux™, OrthoLux™, VertaLux™, SpinaLux™, OrthoLux Lite™, TorsoLux™, Artisan™.
 - Included in every Extender kit is a drill template, instruction kit and Allen wrench for quick attachment to a body jacket.
 - One-size-fits-all, with adjustable height.
 - Each of the 6 sternal and/or pectoral designs (shown here) is available with any of the options below:
 - Fixed or adjustable extension
 - Black, white, or no shoulder straps
 - With soft covers, or only white foam padding.
- Note: All combinations of styles and options are available (see chart below).*
- The "Mickey" and "Minnie," (described below) are available accessories.



EXT PV
Black Straps
Soft Covers




EXT PS
White Straps



EXT STH
White Straps



EXT PR

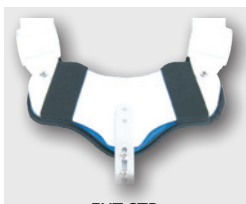

EXT STP
White Straps
Soft Cover



EXT STF

P/N	Style	P/N	Style
07-1031	STF	07-1034	PS
07-1032	STP	07-1035	PR
07-1033	PV	07-1081	STH

Ordering Guide

Orthosis	Style	Extension	Straps	Padding
EXT	STP	ADJ	BLK	WSC
	PR	FXD	WHT	NSC
	PS		NS	
	PV			
	STH			
	STF			



EXT STP
Adjustable Extension

HOPE Orthopedic's "Mickey" posterior over-the-shoulder harness comes as a complete kit with straps & axilla pads. The "Mickey" attaches by Velcro to the Velcro loop back panel of products such as the OrthoLux™ & OrthoLux Lite™, and to Velcro loop patches affixed to the posterior plastic shells of products such as the VertaLux™, SpinaLux™ & ProLux™. (These Velcro patches are available upon request, at no additional charge.) The straps run from the Mickey, anteriorly over the shoulders, posteriorly back under the axilla, and re-attach to the Mickey. A variation which uses the Extender's straps, the "Minnie" is also available.

"Mickey"

